## Sensi Pro

ITO Glass Smart Scale with Body Composition and Sensitive Sensors

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---

## **OUR STORY**

*Born in 1995 in Istanbul, ttec is a global technology brand offering a wide range of well-designed, quality and innovative products at affordable prices with outstanding customer service to beautify people's lives.*

*ttec meets with its customers in 22.000+ stores in 20+ countries. The number of ttec products reached their customers worldwide exceeds 100 million.*



*ttec is owned by Tesan Iletisim AS, Turkey's value added technology distributor with its innovative vision since 1983.*

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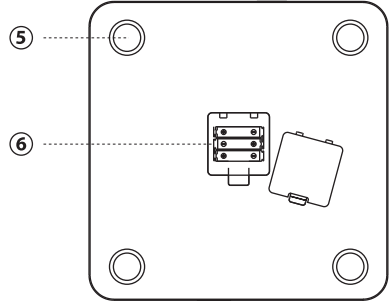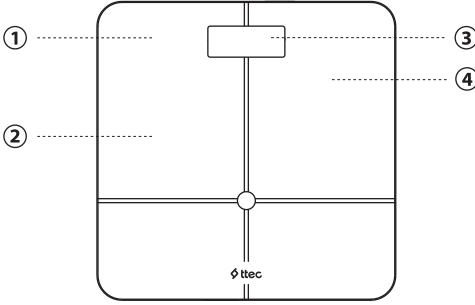## ③ LED Display

LED Gösterge  
LED-Anzeige  
Écran LED  
Pantalla LED  
LED дисплей  
LED дисплей  
Οθόνη LED  
شاشة LED

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## ② Coated Surface with ITO Technology

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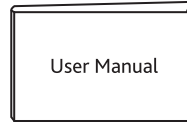
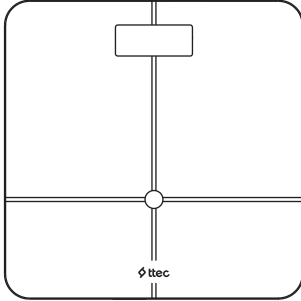
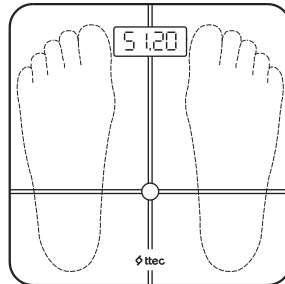
**Carefully step on the platform of the scale. Stay as still as possible until the LED display stabilizes and the weight reading is displayed and make sure that your feet do not stay off the glass surface.**

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### **The scale will be automatically off if without further operation.**

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### **If you don't connect Bluetooth and stand on the scale, the scale will only measure your body weight.**

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### **Please use the scale on the hard & flat surface. Don't use on the carpet or soft surface.**

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## How to connect with Bluetooth for weight, BMI (body mass index), body fat, muscle mass, water, fat-free body weight, bone mass, basal metabolism, protein rate, visceral fat index, body age, body type, body score measurements?

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**• The scale is for family use only and not suitable for professional use. Weight, BMI (body mass index), body fat, muscle mass, water, fat-free body weight, bone mass, basal metabolism, protein rate, visceral fat index, body age, body type, body score measurements are for reference only. You should consult a doctor when you undertake any diet or exercise program.**

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**Caution: Slippery when wet! The scale platform maybe pretty slippery when wet, thus please make sure both the scale platform and your feet are dry before use. Never step on the scale platform with wet foot.**

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Maximum Capacity	180kg/396lb
Auto on and off	Available (The scale would be switched on when weight of 5 kg or more is placed on it.)
Measurement	Weight, BMI (body mass index), body fat, muscle mass, water, fat-free body weight, bone mass, basal metabolism, protein rate, visceral fat index, body age, body type, body score
CPU	High precision strain gauge sensor system and high-performance CPU, make sure of good accuracy.
Tempered glass platform	Available
Low battery and overload indication	Available
Power	3pcs AAA battery with works

**For the first time use, please open the battery cover to insert the battery according to the polarity. If the scale fails to operate, remove the battery completely and reinsert it.**

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**TESAN**  
İLETİŞİM A.Ş.

This symbol means the product must not be discarded as household waste, and should be delivered to an appropriate collection facility for recycling. Proper disposal and recycling helps protect natural resources, human health and the environment. For more information on disposal and recycling of this product, contact your local municipality, disposal service, or the shop where you bought this product.

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